


BONCAST[®]

Synthetic Orthopaedic Casting Tape

Fiberglass and Polyester Substrate

BONCAST[®] is a synthetic orthopaedic casting tape made of knitted fiberglass or polyester substrate impregnated with water-activating polyurethane resin. The fabric is extensible in all directions, which makes it very comfortable in use.

Exposure to any moisture will cause the tape to become rigid. The setting time is approximately 4 minutes and it takes 20 minutes for full weight bearing.

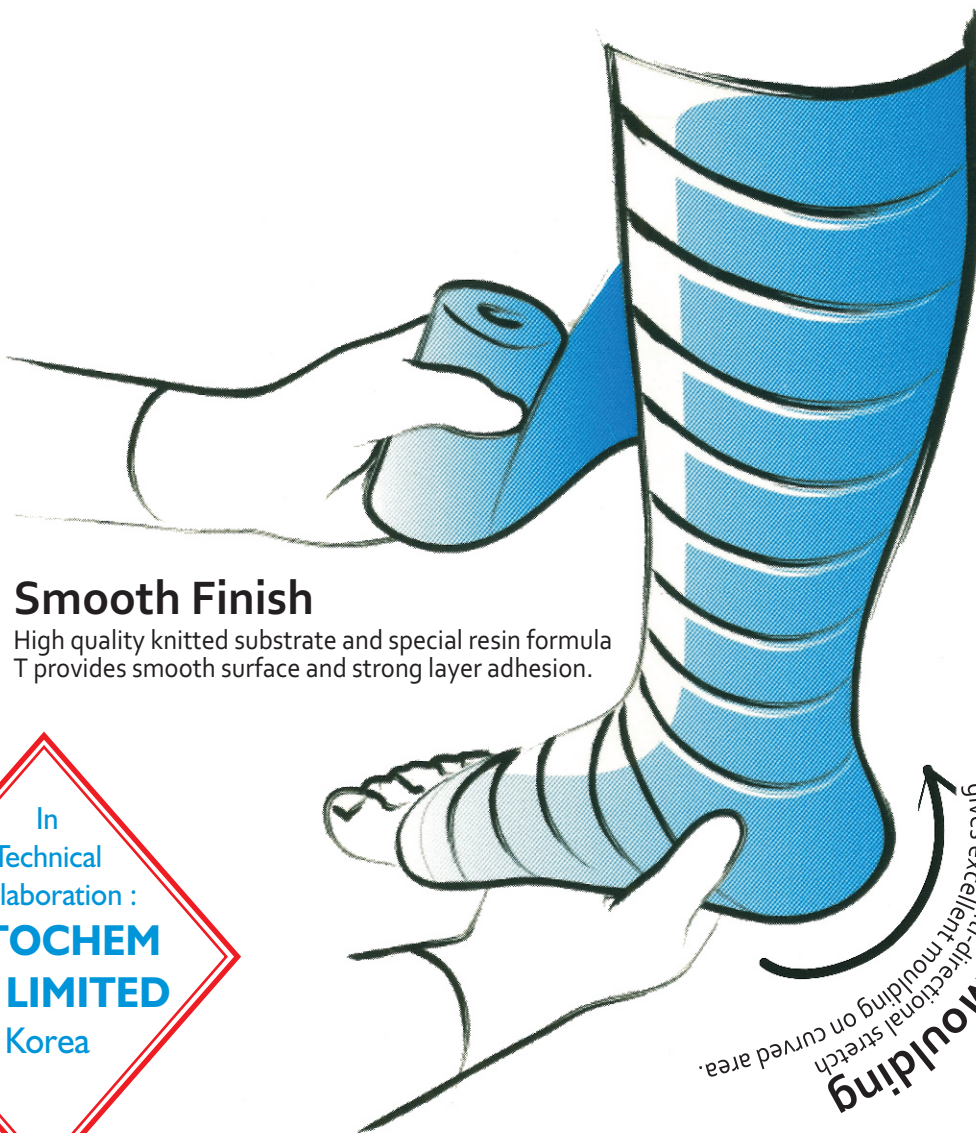

Smooth Finish

High quality knitted substrate and special resin formula T provides smooth surface and strong layer adhesion.

In
Technical
Collaboration :
**ALTOCHEM
CO. LIMITED**
Korea

Excellent Moulding
Superior multi-directional stretch
gives excellent moulding on curved area.



DERMACARE SURGICALS PVT. LTD.,

Plot No. – 120, Village Mahim, Taluka Palghar, Dist. Thane – 401 404.

Email : info@dermacare.co.in Visit us at : www.dermacare.co.in

Customer Care No : +91 9323012321



BONCAST[®] IS...

Light weight

Excellent Moulding

Strong

Superb Adhesion

Air Permeable

Smooth Unwinding

Moisture Resistant


Simply to use

X-ray Transparent



Full range

No End-edge problem

It provides No End-edge folding problem,
which remains doctors time consuming and in convenience job of unfolding its end.

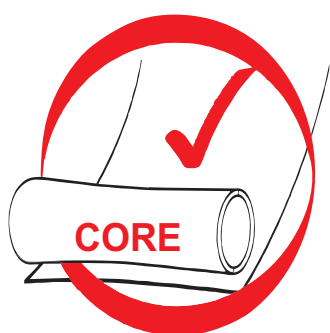

BONCAST[®]

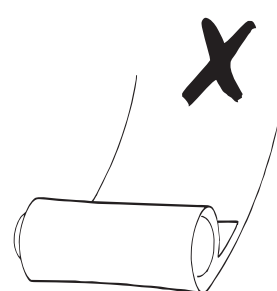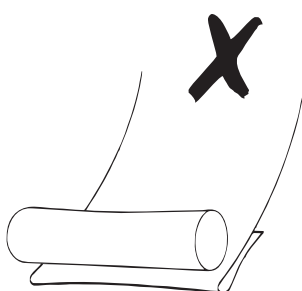

Others

CODE NO.	SIZE	PACKAGING	SUBSTRATE	COLOUR
DF 7002	2in x 4yds(5.0cm x 3.6m)	1 Roll / Pouch 10 Rolls / Inbox 100 Rolls / Outbox	FIBERGLASS / POLYESTER	Blue, Green and Pink, (Other Colours on request)
DF 7003	3in x 4yds(7.5cm x 3.6m)			
DF 7004	4in x 4yds(10.0cm x 3.6m)			
DF 7005	5in x 4yds(12.5cm x 3.6m)			